

# ANESTHESIA CARTS

Organization and work surface are important features on an anesthesia cart and the Waterloo aluminum anesthesia cart offers both. With the U-Divide drawer divider system included in every drawer a Waterloo aluminum cart is organized right out of the box. Fine tune your organization with a vast variety of custom trays and dividers that offer more specific organizational options.

With two pull out side shelves, work surface on a Waterloo aluminum cart is virtually double that of competitive carts! Couple that with a variety of above cart tilt bin and shelving options and you have all the flexibility needed to organize equipment and supplies that need to be readily available in the OR.



UTGSA-333369-DKB

Dual slide out work surfaces operate independently and virtually double your work surface!

Included plastic top with cutaway front is designed to be "Wrist friendly" as you work on the cart top.

Large stylish/contemporary drawer handles offer a more updated look and enhanced functionality


Adjustable metal drawer dividers (One long and One short) included in all aluminum cart drawers.

Advanced nylon reinforced polycarbonate bumper provides a wider wheel base and greater cart stability and protection. Single piece construction resists scuffs and scratches.

5" Premium "Easy roll" sealed bearing casters with precision bearing swivels. Thread guards minimize sutures and other materials from collecting around axes and keep carts rolling smoothly.

RECOMMENDED

**TILT BINS**  
 Convenient tip out tilt bins available in 3Bin, 4Bin, 5Bin and 6Bin sizes. Requires Base rail system sold separately  
 Part#s TB-3, TB-4, TB-5, TB-6  
 SEE PAGE 321

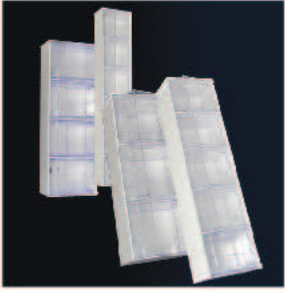

**BASE RAIL SYSTEM**  
 Horizontal rail system for mounting tilt-bins, shelves and other accessories. Part#s BRS-2 or BRS-Jr (Med Jr Carts) SEE PAGE 32.



**WASTE BASKET**  
 Large 20Quart plastic waste basket. Part#s BSK-2 (20Qt) or BSK-1 (6Qt) SEE PAGE 31.



FEATURE

ELITE ALUMINUM

Unicart Construction

Ball Bearing Slides

Easy-roll Casters

Includes Drawer Divider Set

Universal Accessory Mounting System

Dual Pullout Work Shelves

Lightweight Construction



**CLEAR ORGANIZER**

Clear acrylic organizer to keep small items at hand.  
Part# ORC-1S  
SEE PAGE 32.



**TAPE DISPENSER**

Light gray metal tape dispenser can be mounted directly to the top of the cart via included suction cups or to a base rail system using the included accessory hooks. Holds approximately 10 rolls of tape.  
Part# TD-1  
SEE PAGE 32.



**WIRE STORAGE BASKET**

Extra Large Wire Storage Basket mounts to accessory railing and measures 15.5"W x 9.5"D x 5.5"H.  
Part# WSB-1  
SEE PAGE 32.

**COMPARISON**

**TRADITIONAL STEEL**

- Unicart Construction
- Ball Bearing Slides
- Easy-roll Casters
- Includes Drawer Divider Set
- Universal Accessory Mounting System
- Dual Pullout Work Shelves
- Lightweight Construction

**ANESTHESIA CARTS**

Economical steel anesthesia carts offer a robust work platform that can be configured with accessories to facilitate your work in the OR.

Shares many of the same great features as our aluminum cart but at a more economical initial price.

Collapsible side shelves are available for added work surface and a variety of drawer trays and dividers are offered to further organize your cart.

Choose from a variety of trays and dividers to organize critical medications and supplies.



If you are considering a steel cart and are looking to organize it with plastic trays/dividers as well as a side shelf or two, you may wish to consider our aluminum cart as the included dividers and side shelves can offset the added cost of upgrading to aluminum.



Smooth, easy-to-clean and removable work surface.

Pre-drilled holes to mount accessories

Advanced nylon reinforced polycarbonate bumper provides a wider wheel base and greater cart stability and protection. Single piece construction resists scuffs and scratches.

5" Premium "Easy roll" sealed bearing casters with precision bearing swivels. Thread guards minimize sutures and other materials from collecting around axles and keep carts rolling smoothly.

UTGSU-433369-DKB

## ANESTHESIA CARTS

Anesthesia carts come in all shapes and sizes, and with various security requirements. At Waterloo we acknowledge that fact by manufacturing a wide variety of cart sizes and locking options designed to fit your specific needs. Simply choose your desired drawer combination of 3", 6", 9" or 12" drawers and we can help you pick the correct cart height to accommodate your configuration.

Then simply select your cart exterior color and your drawer front colors. If you'd like to color code drawers to make them easy to identify quickly (example: White airway drawer), we are happy to accommodate that request at no extra charge! Just tell us what drawer front and what color and we can personalize your carts for you. Another advantage of working with a true manufacturer of medical carts!

Waterloo aluminum anesthesia carts are available in Extra-Tall, Tall, Medium and Short cart sizes (See page 5 for sizing chart).

Waterloo anesthesia carts are available with a variety of locking options. From a simple key lock to a WiFi capable proximity card reader, we can provide the right level of security for your operating rooms.



Further enhance cart security with an independently locking "controlled substance" drawer (optional) to keep drugs secure. One swipe of an approved employee badge gains access to all drawers with the exception of the top drug drawer. A secondary key-code allows access to the controlled substance drug drawer keeping important drugs secure.






## ANESTHESIA CARTS

Steel carts come in all of the same cart heights as our versatile aluminum carts.

As with our aluminum carts, choose your desired drawer combination of 3", 6", 9", or 12" drawers. From there we can determine the appropriate cart height for your configuration. For available cart sizes see sizing chart on page 5.

Next, select your cart exterior color and preferred drawer front colors. Should you like to color-code any drawers, simply let us know which drawer and what color (from one of our 16 available colors) and we are happy to accommodate your request. We never charge for this option!

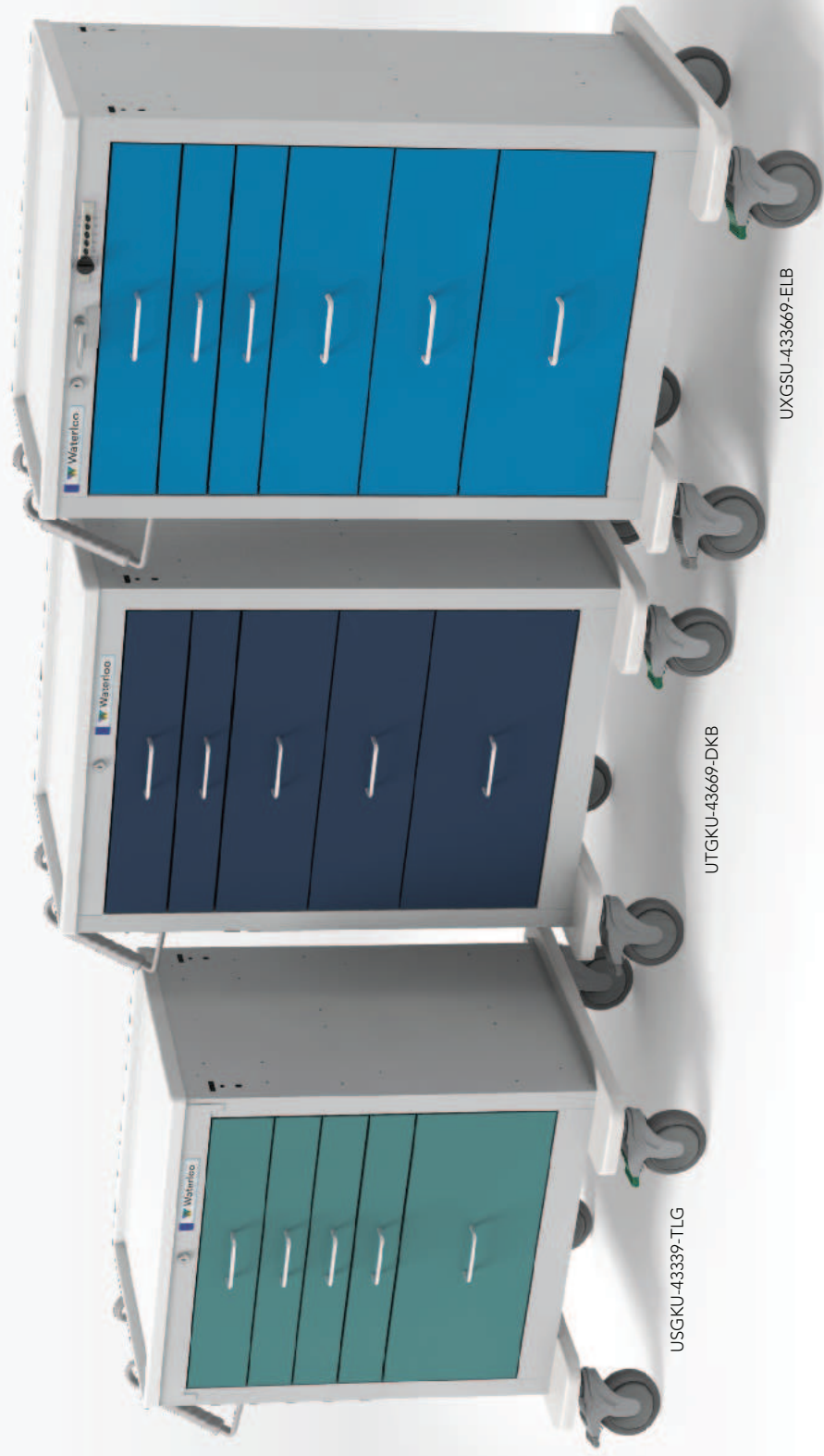

USGKU-43339-TLG

UTGKU-43669-DKB

UXGSU-433669-ELB