NewCartago Armchair





Features

The armchairs are specially design to be used in hospitals, clinics, nursing homes and also in home care.

Their easy to use makes it robust and resistant for extensive uses.

The armchair is made of a tubular steel structure with polyester coating, with upholstered seat, backrest and footrest sections with multiple colour options.

The armchair is articulated in two sections by levers located in both sides, driven by gas spring to raise and low the backrest and footrest.

It also has foldable armrest to allow side access.

Depending on the intended user several configurations are available based on the quantity of wheels it has or the accessories can be placed in it.

Options:

For patient use, it can integrate 2 or 4 wheels, pushing and carrying handle, table, tank holder, IV pole and footrest.

For companion use, no wheels or other element will be necessary.













Seat height











NewCartago Armchair

pardo

Technical information


Total length (seated position) (A)	850 mm.
Total length (tilted position) (B)	1610 mm.
Total width with armrests (C)	605 mm.
Total width (D)	1235 mm.
Seat height (E)	490 mm.
Seat width (F)	520 mm.
Backrest tilting angle (relative to the seat back) (G)	118°
Extremities tilting angle (relative to the seat back) (H)	88°
Maximum backrest tilting angle (relative to the seat back) (I)	155°
Maximum extremities tilting angle (relative to the seat back) (J)	9°
Armchair weight	24 kg.
Maximum patient weight (MPW)	135 kg.






IV Pole


Upholstery

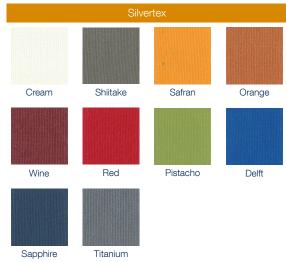



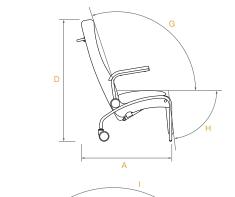



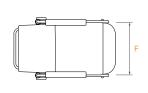



- 4 legs.
- Without handle.
- Without accessories.



versión with two wheels

- 2 wheels.
- 2 wheels Ø 100 mm.
- With handle.
- Optional accessories.



tient versión with four wheels

- 4 wheels Ø 100 mm
- With handle.
- With footrest.
- Optional accessories.